

YSI 9800 reagent instructions

Fluoride (YAP179L)

Photometric test for Fluoride in drinking water, natural water, and treated water

Color change

Red to yellow



Range

0 - 1.5 mg/L F

01



Fill to the 10 mL line with sample.

02



Use this to "blank" the instrument. Place in the cell holder and press **Blank**.

03




Prepare the sample by adding 4 drops of the Fluoride liquid reagent and stir to mix.

04



Add a Fluoride No 2 tablet, crush, and stir.

05

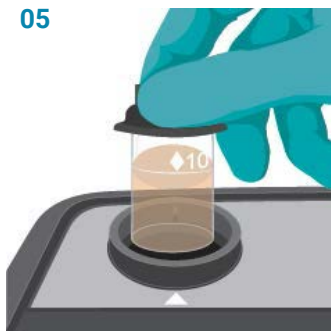

Place in the cell holder and press **Measure**.

06



Allow the timer to count down 5 minutes. The concentration of Fluoride will then be displayed.